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TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY. — Nebezpečenstvo udusenía.

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Symbol  znamená UPOZORNENIE.

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Sonic toothbrush**eta 7707****SONETIC****INSTRUCTIONS FOR USE**

Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.

**I. GENERAL SAFETY INSTRUCTIONS**

- Consider the instructions for use as a part of the appliance and pass them on to any other user of the appliance.
- Check whether the information on the name plate corresponds to the voltage in your electrical socket. The power supply must only be connected to the corresponding AC mains adapter or backup power supply corresponding to the relevant technical parameters.
- This product is not intended for children under 3. Children at the age of 3 to 14 and persons with limited physical, sensory or mental capabilities may only use the toothbrush if they are supervised or have been instructed of the use of the appliance in a safe manner and understand possible risks. Children must not play with the appliance. Children must not perform cleaning and maintenance intended to be performed by users, unless they are older than 8 years and are under supervision. Children under 8 years of age must be kept out of reach of the appliance and its charging base.
- The appliance has no charging source of its own, it is only designed for charging via USB (an adapter with USB output with output voltage of 5.0 V, and output current of at least 600 mA). Before making any connection, please pay attention to instructions of the manufacturer of the respective power source.
- The connected charging base must not be used at the places where it could fall in bathtub, sink or swimming pool. If it falls into water yet, do not take it out! First disconnect the appliance from the electrical socket, and then take it out from water. In such cases bring the application to the repair shop to verify its safety and correct functions.

- Never use the appliance if the charging base or power cord is damaged, if it operates incorrectly, fell on ground or is damaged, or if the charging base fell into water. In such cases bring the appliance to the specialised repair shop to verify its safety and correct functionality.
- To power the appliance, only apply safe small voltage of 5V.
- The appliance is only intended for individual use at home for personal oral cavity care! It is not intended for commercial use, for more than one patient in dentist's consulting rooms or in any other healthcare facilities!
- The charging base provided with appliance is only designated for this appliance. Do not use it for other purposes. Only use the provided charging base for this appliance, do not use other types of charging base for recharging (for example from other appliance).
- Do not use the application in exterior or where aerosol products are used (such as **sprays**), or where oxygen is handled.
- **Do not use the charging base during bathing.**
- Guarantee of 24 months does not apply to decrease in the accumulator's effective capacity caused by its use or age. Effective capacity gradually decreases depending on the method of use of the accumulator.
- Do not insert the lead of the charging base into the USB power supply or remove it with wet hands!
- **Never immerse the charging base in water or other liquids (not even partially) and do not wash it under running water!**
- Do not insert the appliance or its accessories in any body cavities, except for oral cavity.
- The appliance is not designed to care for animal teeth.
- Only use the charging base at horizontal position at places where overturning is not imminent, and in a sufficient distance from heat sources (such as **stove, fireplace, or heat radiators**).
- Clean the accessories after each use. Failure to adhere to this instruction might cause worsening your state of health, because without disinfection, spread of mildew and bacteria is supported.
- Charge the accumulator at the normal room temperature. Do not expose the accumulator to temperatures higher than 50 °C to prevent causing damage to the accumulator.
- The accumulator must not be removed when the appliance is normally used. Only remove the accumulator if the appliance is to be discarded (see article **VI. ECOLOGY**).
- **Do not leave the appliance unattended!**
- Protect the appliance from dust, chemicals, toiletries, liquid cosmetics, large fluctuations in temperature, and direct sunshine.


- Do not attempt to repair the appliance by yourselves.

SPECIAL SAFETY INSTRUCTIONS

- If you have unpleasant feeling in your mouth or if you feel pain during the cleaning, stop using the appliance and contact your dentist. Contact the dentist also in the case that you have any doubt about use of the appliance in your individual case.
- If your toothpaste contains peroxide, baking soda or other bicarbonate (frequently used in whitening toothpastes), thoroughly clean the brush and its handle after use to prevent damaging the appliance.
- Do not use the appliance in care for children, incapacitated persons or persons who are unable to express their feelings.
- **If you have undergone tooth or gum operation recently, suffered long-term problems with your teeth or were subject to extensive medical interventions in the last 2 months, or if you suspect problems of dental character (with fillings, crowns, etc.), consult your dentist before using the appliance.**
- **If you experience any health troubles (for example considerable bleeding, or if any bleeding persists even after 2 weeks of use of the brush), stop using the appliance and contact your physician for consultation.**
- The appliance complies with the safety standards for electromagnetic radiation. However if you use any of the following devices: cardiostimulator, electrocardiograph or other medical device implanted in your body or used as a support, consult your physician or the implant manufacturer before use.
- For hygienic reasons do not share the toothbrush with another person or more persons. Infections could spread or oral cavity inflammation could be initiated.
- The charging base power cord must not be damaged by sharp or hot items or open flame, and they must not be immersed in water or bent over sharp edges. They must not freely hang over the desktop edge within reach of children.
- Regularly check the condition of the charging base power cord.
- Only use the appliance with the original accessories designated for the respective type. Do not use damaged/deformed adapters.
- Only use this appliance including its accessories for the designated purpose as described in this instruction manual. The appliance is only designed for cleaning teeth, gums and tongue. Never use the appliance for any other purpose.
- **WARNING:** There is a risk of injury in case of incorrect use (not in accordance with the manual) of the appliance.
- All text in other languages and pictures shown on the packaging, or the product itself, are translated and explained at the end of this language mutation.
- The producer does not undertake any responsibility for damage caused by improper use of the appliance (e.g. **injury etc.**) and its guarantee does not cover the appliance in the case of failure to comply with the safety warning above.

Wireless connectivity and compatibility

- ETA can in no way be held liable for any accidental, indirect or subsequent damage, data loss or any other loss incurred as a result of information leakage during communication through BLUETOOTH, Wifi technology.
- 100% successful wireless communication cannot be guaranteed with all Smart devices. Devices have to be compatible with BLUETOOTH specifications set forth by the company Bluetooth SIG, Inc., and their compliance with the specifications has to be certified. However, given the variety of devices available in the market, even if a device is in compliance with the specifications, there might be cases when the character or specifications of the given Smart device do not allow connectivity, or hinder connectivity in any way. This is not a reason to make a warranty claim.
- Radio waves transmitted by devices may interfere with operation of medical equipment. Such interference may cause a failure. Therefore, do not use the device in the following places:
 - in hospitals, trains, planes, petrol stations and in areas where flammable gases may be present.
 - in the vicinity of automatic doors or fire detectors.

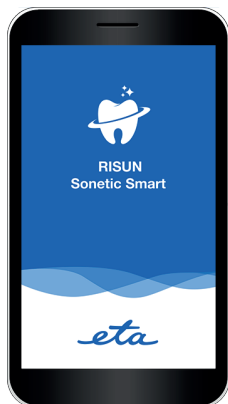The application is owned and operated by a third party (hereinafter the “Subject”) different from HP TRONIC Zlín, Ltd. (hereinafter the „Seller“). No personal information is transmitted between the Subject and the Seller. When registering for the application, the customer transfers the personal data directly to the Subject, which becomes the personal data controller. The Application and the associated privacy policy are governed by the Rules of the Subject.

Remark

The application requires an active internet connection.

II. GENERAL INFORMATION

You have bought a device, at a first glance looking and acting as an ordinary toothbrush. However, if connected via bluetooth to any Smart device with Android or iOS operating system with RISUN App installed, the toothbrush becomes a multifunctional modern assistant for monitoring and better evaluation of your dental care.

App features

- better overview of dental hygiene.
- advise on which hygiene cleaning plan to choose according to the condition of your teeth.
- brush cleaning and usage statistics.
- possibility of individual setting of intensity and cleaning time.
- recommended cleaning techniques and procedures.
- brief evaluation of the cleaning technique.
- practical information on oral hygiene.
- possibility of history sharing.
- cleaning reminder function (so-called alarm clock).
- Designed for Android and iOS (can be downloaded for free from Apple Store / Google Play),

First steps

Before first use, download the **RISUN** application to your Smart device for free. You can find the application in the Apple Store or Google Play database, or to speed up the process, you can use the following QR codes:

III. DESCRIPTION OF FEATURES

A – Brush body

- A1 – ON/OFF switch + mode switch
- A2 – wireless connection indicator ●
- A3 – battery indicator light
- A4 – user mode indicator
- A5 – normal mode
- A6 – whitening mode
- A7 – intensive cleaning mode
- A8 – massage mode

B – Brush (3 pieces included)

- B1 – brush head
- B2 – coloured ring to distinguish users

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IV. INSTRUCTIONS FOR USE

Remove all packaging materials and take out the appliance and any accessories. Parts which come in contact with the oral cavity must be rinsed in warm drinking water before first use.

CHARGING THE BRUSH

This appliance is equipped with a contactless (induction) charging system which works on the principle of energy transfer through an electromagnetic field. No direct connection with the power cable is required which makes charging more convenient and thanks to the enhanced insulation (IPX7) also safer. First, connect the charging base to a suitable power source using a USB power cable. Then place the toothbrush body (**A**) on the cylindrical protrusion in the charging base (**C**). The charging procedure is signalled by a flashing (**A3**) LED indicator. When the battery is fully charged, the (**A3**) LED indicator is lit. For a full charge, the battery needs to be charged for approximately 12 hours. Depending on the cleaning mode used, a fully charged battery will provide approximately 4-5 hours of operation.

Note: Although the brush should not be charged before it is fully discharged, we recommend charging it at regular intervals (e.g. after 2-3 weeks of use) in order to ensure optimum battery life.

USE

Before use, slide the brush (**B**) onto the brush body (**A**) to connect them (See Fig. 1). A small gap will remain between the brush body (**A**) and the coloured ring (**B2**) which is not a defect. To remove the brush, proceed the other way round. Do not switch on the appliance without the brush attached!

Moisten the brush and apply a small amount of toothpaste on it. Put the brush to the teeth at a slight angle towards the gums (see Fig. 3) and switch on any mode by pressing the mode on/switch button (**A1**). Move the brush slowly along the teeth back and forth, applying a slight pressure, so as not to miss any tooth or area in the mouth. Clean each tooth for about 2-3 seconds.

The operating time of the brush may be different for each mode. It can also be changed in the application. After this time, the brush switches off. If you want to continue cleaning, simply turn the brush on again.

According to general recommendations, it should take approximately 30 seconds to clean each part of the mouth (quadrant) (see Fig. 3); however, this time is individual and depends on the settings and mode used. A short interruption of brush operation will alert you to the expiration of a quarter time and a change in the quadrant.

Caution

- During operation, the (A3) LED indicator is lit, when the battery of the accumulator is low and it needs to be recharged, the LED indicator starts flashing.
- When charged for the first time, or if the appliance was not in use for more than 3 months, it needs to be charged for 12 hours.
- **Do not exceed the recommended time of charging during regular charging!**
- The charging base also serves as a stand on which you can place the toothbrush after use. However, do not connect it to an external source unless you want to charge the toothbrush. Do not leave the charging base always connected to an external source!
- The appliance warms up during charging, but this is normal and not a defect.
- The operating time is only approximate. The time depends on the condition of the battery, ambient temperature, age, method of charging, etc.
- To power the appliance, use only a safe low voltage source of 5 V DC in accordance with the applicable standards (a mobile phone charger with a USB connector in double insulation is suitable), which is capable of giving a current of at least 600 mA or more. (Note: the specified current load is minimal, it can be higher).
- Do not use USB ports on other devices (e.g. computer, laptop, tablet, mobile phone, etc.) to power the appliance. These devices do not have to be designed for current consumption of the appliance and could be damaged.

Notes

- If this is your first sonic toothbrush and you have not had any experience with this cleaning technology, it might take you a while to get used to the specific vibration (tingling). However, most users say that they got used to it quickly after a few uses and cleaning became perfectly natural for them.
- Do not apply excessive pressure on your teeth.
- You can also use this brush when you have fixed braces, facets, crowns, etc. However, please bear in mind that if the brush is used with braces, the brush bristles are subject to more wear.
- We recommend that you replace the brush (B) with a new one after 3 months of use. However, if the brush is showing signs of significant wear or damage.

Mode on/switch button control A1

The first press of the button switches the brush to the last mode used. If you press the button again within 3 seconds, the brush will switch to the next mode. The last press turns the brush off. The brush also turns off if you press the button (**A1**) even in the set mode 3 seconds after its activation.

Cleaning modes:

- 1) mode **A5** – standard intensive mode for routine cleaning of healthy teeth,
 - 2) mode **A6** – cleaning mode adapted for better teeth whitening,
 - 3) mode **A7** – mode for polishing and more intensive cleaning,
 - 4) mode **A8** – pulsation mode for gentle massaging of the gums,
 - 5) **A4** mode – individual mode for custom cleaning settings in the application.
- Activation of the mode selected is indicated by the indicator coming on.

APP

The toothbrush can be used with the Smart application RISUN. After finding and installing the App, register an account, agree with necessary conditions, and enable necessary permissions for the App to operate properly.

Remarks

- If you are using VPN on your phone, please turn it off to sign in and register. This will avoid any problems with logging into the server. Some other features of the application may be affected by using of VPN. If you have a problem with them, disable it to work with the application.

Pairing and use

To pair the brush with the App for the first time, click on pairing at the top of the screen or directly in the settings in the My brush section. Then follow the instructions provided directly in the application. The App automatically finds the brush nearby. After selection, the wireless connection indicator **A2** will start flashing on the brush for a few moments during which it is necessary to press the button **A1** to successfully complete the pairing.

Then the application will offer you the possibility to adjust the settings and the recommended cleaning mode according to the condition of your teeth. You can skip this option and use it at any time later in the relevant part of the App.

The App will display and record the way the brush is used, at the same time you can adjust some parameters here. When using, turn on the application first, then the toothbrush. Alternatively, you can use the brush without the App - after activating it, the data is automatically synchronized and the data from the paired brush is uploaded to the App.

If the App and the brush are used by multiple household members, each member should use their own account with its own settings and history. At the same time, no more than one application can be activated at a time. However, when used by more than one person, there may be various inaccuracies in the data synchronization (for example, if you do not use the brush with the App and someone else has paired it with your account in the meantime, your cleaning history may not upload correctly or your data may be uploaded to another account).

If the devices do not communicate

From time to time, pairing with the Smart device might fail. If this is the case, try one of the following steps, or a combination thereof:

- 1) Place the devices as close as possible to each other (approx. 10 cm)
- 2) Make sure that the system is not disturbed by another el. appliance! (Eg. remote control TV, etc.)!
- 3) Restart the appliance.
- 4) Restart the **RISUN** application, or switch it off and on in the system.
- 5) Restart the entire Smart device.
- 6) Try the application with another device.

Note:

- Make sure you have the bluetooth function activated on your Smart device with permanent visibility to other devices, internet turned on and the App has all necessary permissions enabled (e.g. localization function if required).
- If you are unable to pair the App even after trying all the advice given in the manual and the link and the brush does not respond to other compatible Smart devices, only then make a complaint about the brush.
- Working with the application is intuitive and it won't take you long to control it. Experiment, discover and enjoy the possibilities.

V. CLEANING AND MAINTENANCE

Before servicing the charging base (C), first disconnect it from the external power supply. For maximum hygiene, clean the brush (B) after each use! Rinse the brush head (B1) and the brush body (A) under running water and wipe the brush body with a damp cloth. Use the damp cloth to wipe the charging base and its cable as well. Never immerse the brush body (A) or the charging base (C) in water! Do not use a dishwasher for cleaning! Allow the brush body (A), the charging base (C) and the cable to dry properly before re-use.

Store the appliance and all the accessories in a dry, dust-free and safe place, out of reach of children and incapacitated persons.

VI. ENVIRONMENTAL PROTECTION

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal of the product, hand it over at special collection places where it will be taken over free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. This appliance is equipped with a Li-ion battery with long life. For environmental protection it is necessary to remove the battery from the appliance after its service life and the battery and the appliance must be disposed of safely using special collection networks. Ask for more details at the local authority or at the local waste collection site.

Removing the battery

Remove the battery from the appliance only when completely discharged.

The procedure for removing the battery is shown in Fig. 4:

1. Using a screw-driver, remove the screw cover and loosen the screw.
2. Hold the brush body (A) upside down, press the charging part (see arrows in the figure) and remove the 2 hooks which hold the inner components of the brush body (A).
3. Remove the battery from the compartment using a screw-driver.

VII. TECHNICAL DATA

Voltage (V) appliance	shown on the type label of the
Weight (kg) approximately	0.12
Appliance/Adapter Protection Class	III./II.
Size, (mm):	250 x 33 x 33
Frequency range:	2,4 - 2,483 GHz
Output power (max):	0 dBm
Bluetooth version:	4.2
Compatible operating systems:	Android 4.3; iPhone4S a iOS 8.0 (and higher)
Maximum communication reach:	10 m (open space without any obstacles and interferences)

Hereby, ETA a.s. declares that the radio equipment type ETA7707 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the ETA www address.

The manufacturer reserves the right to change the technical specifications and accessories for the respective models. Expert service must do the large maintenance or to intervene into inside parts of the appliance! Non-observance of the instructions of the manufacturer cancels the right for guarantee repair!

The manufacturer reserves the right to change/update the RISUN application without prior notice. Due to the character of such changes/updates, some information or procedures in this manual might become obsolete or might be missing.

WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:

HOUSEHOLD USE ONLY; DO NOT IMMERSE IN WATER OR OTHER LIQUIDS;

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

The  symbol indicates a WARNING.

- Please read the instructions manual

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© DATE 6/8/2020