05584

DE

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ab Seite 3

EN

Operating Instructions

from page 11

FR

Mode d'emploi

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NL

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vanaf pagina 27

Inhalt

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
Contents

Intended Use	12
Safety Notices	12
Items Supplied	15
Device Overview	15
Use	16
Changing the Battery of the Remote Control	16
Cleaning and Storage	17
Troubleshooting	17
Technical Data	18
Disposal	18

IMPORTANT NOTICES! PLEASE KEEP FOR REFERENCE!

Explanation of the Symbols

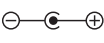-  Danger symbols: These symbols indicate risks of injury (e.g. from electricity or fire). Read through the associated safety notices carefully and follow them.
-  Read through the associated safety notices carefully and follow them.
-  Read through the associated safety notices carefully and follow them.
-  Only use indoors!
-  Polarity of the coaxial connector (mains adapter)
-  Symbol for direct current
-  Protection class II
-  Safe electrical isolation
-  Switching power supply
-  Supplementary information

Read operating instructions before use!

Explanation of the Signal Words

- DANGER** warns of serious injuries and danger to life
- WARNING** warns of *possible* serious injuries and danger to life
- CAUTION** warns of slight to moderate injuries
- NOTICE** warns of material damage

Dear Customer,

We are delighted that you have chosen this **massage device**.

Before using the device for the first time, please read the operating instructions carefully and keep them for future reference. If the device is given to someone else, it should always be accompanied by these operating instructions. They are an integral part of the device.

The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

If you have any questions about the device and about spare parts / accessories, contact the customer service department via our website:

www.dspro.de/kundenservice

Intended Use

- The device is intended for massaging and stimulating your feet and leg muscles and is used for relaxation. It is not a replacement for medically required treatment.
- The device is only suitable for use when sitting or lying down.
- Do not use the device before going to bed. The massage has a stimulating effect, which may lead to problems getting to sleep.
- The device is not designed to be operated with an external timer or a separate remote control system.
- The device is for personal use only and is not intended for commercial applications.
- The device should only be used as described in the operating instructions. Any other use is deemed to be improper.
- The warranty does not cover faults caused by incorrect handling, damage or attempted repairs. The same applies to normal wear and tear.

Safety Notices

- This device should not be used by people (including children) with impaired sensory or mental faculties or who do not have adequate experience and / or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.
- **Children** must not play with the device.
- Cleaning and **user maintenance** may not be performed by **children**, unless they are supervised.
- Do not use the mains adapter with another device and never use the device with another mains adapter.

- Check the device for damage each time before you use it. The device may only be used if the device and connecting cable have no visible signs of damage! If the device or the connecting cable are damaged, they must be replaced by the manufacturer, customer service department or a specialist workshop in order to avoid any hazards.

Dangers for the Following Groups of People

- The following people may use the device **only after consultation with their doctor**:
 - people who wear a pacemaker or similar medical implants;
 - people with heart disease or who suffer from epilepsy;
 - people who have back complaints or spinal disorders;
 - people who suffer from osteomyelitis, bone inflammation (in particular if the spinal column is damaged) or from osteoporosis;
 - women who are pregnant or have just given birth;
 - people who suffer from circulation problems or thrombosis;
 - people who have undergone an operation for which the treatment is not yet complete;
 - people with physical abnormalities;
 - people who have a malignant tumour.
- If you are unsure whether the device may be used, consult a responsible doctor.
- Do not use the device near to or on injuries, swelling, skin irritations, infected skin, warts, fungal infections, abnormal changes in the skin or body parts with no feeling.

DANGER – Danger of Electric Shock

- Use and store the device only in closed rooms. Do not operate the device in rooms with high humidity.
- Never immerse the device in water or other liquids! Ensure that the device and the connecting cable cannot fall into water or become wet.
- Should the device fall into water, switch off the power supply immediately. Do not attempt to pull the device out of the water while it is still connected to the mains power!
- Never touch the device and the connecting cable with wet hands when these components are connected to the mains power.
- Do not use the device if it has malfunctioned or been dropped in water. Have it checked in a specialist workshop before using it again.
- The device must always be disconnected from the mains power when left unattended and before cleaning it.

**WARNING – Danger of Injury**

- **Danger of suffocation!** Keep children and animals away from the packaging material.
- Children younger than 8 years of age and animals must be kept away from the device and the connecting cable.
- **Danger of burns from battery acid!**
 - Batteries can be deadly if swallowed. The battery and remote control should therefore be kept out of the reach of children and animals. If a battery is swallowed, immediate medical attention must be sought.
 - If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. Wear protective gloves if necessary. If you touch battery acid, rinse the affected areas with abundant clear water straight away and seek medical assistance immediately.
- **Danger of tripping!** Make sure that the connecting cable does not present a trip hazard.

**WARNING – Danger of Fire**

- Do not use the device in rooms containing easily ignitable or explosive substances.
- Do not cover the device during operation in order to prevent it from catching fire.
- Unwind the connecting cable fully before connecting to the mains power.

NOTICE – Risk of Damage to Material and Property

- Only connect the mains adapter to a properly installed plug socket. The plug socket must also be readily accessible after connection so that the connection to the mains can quickly be isolated. The mains voltage must match the rating plate of the device and the mains adapter.
- The device must be switched off when the mains adapter is pulled out or inserted into the plug socket.
- Make sure that the connecting cable is not squashed, bent or laid over sharp edges and does not come into contact with hot surfaces.
- Switch off the device and pull the mains adapter out of the plug socket if an error occurs during operation or before a thunderstorm.
- When the mains adapter is removed from the plug socket, always pull on the mains adapter and never the mains cable.
- Do not pull or carry the device by the connecting cable. Never wind the mains cable around the device or mains adapter.
- Use only the battery type specified in the technical data.
- Remove the battery when it is dead or if you will not be using the device for a prolonged period. This will prevent any damage caused by leaks.
- Remove a leaking battery immediately.
- Batteries must not be charged or reactivated in any way, disassembled, thrown into a fire, immersed in liquids or short-circuited.

- Use only original accessories from the manufacturer in order to guarantee that there is no interference that may prevent the device from working and to avoid possible damage.
- Never place heavy objects on the device or the connecting cable.
- Protect the device from heat sources, naked flames, sub-zero temperatures, persistent moisture, wet conditions and impacts.


Items Supplied

- 1 x massage device
- 1 x mains adapter
- 1 x remote control incl. 1 x button cell battery (3 V, CR2025; already inserted)
- 1 x operating instructions

Check the items supplied for completeness and the components for transport damage. If you find any damage, do not use the device but contact our customer service department.

Remove any possible films, stickers or transport protection from the device. **Never remove the rating plate and any possible warnings!**

Device Overview

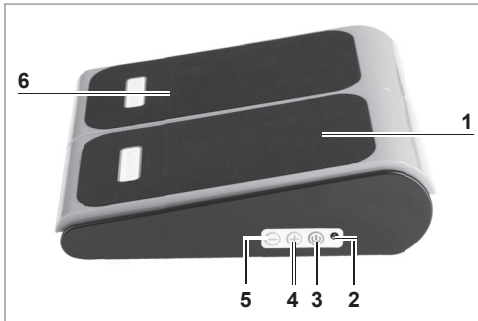- 1 Vibrating surface on right
- 2 Connecting socket for the mains adapter
- 3 Button 
- 4 Button +
- 5 Button -
- 6 Vibrating surface on left

Use

PLEASE NOTE!

- No pressure should be exerted on the device as it may be damaged. Never stand on the device.
- **CAUTION – Danger of Crushing!** Never place your hands and feet between the frame and vibrating surfaces while the device is in operation.
- Do not use the device for longer than 15 minutes without stopping.
- Do not fall asleep while you are using the device!
- A tingling sensation / itching during or after the massage is normal due to the circulation-boosting effect.

1. Pull the insulation strip out of the battery compartment of the remote control to activate the battery.
2. Place the massage device on the floor. The surface should be firm, dry and flat. **Please Note:** Protect any delicate surfaces (laminated, parquet floor, etc.) from scratches etc. with a protective mat.
3. Connect the barrel connector of the mains adapter to the connecting socket (2) on the device and insert the mains adapter into an easily accessible plug socket.
4. The device can be used either while sitting or lying down:


Sitting: Sit down on a chair, armchair or similar place in front of the device. Place your feet down on the vibrating surfaces (1, 6).

Lying down: Lie down in front of the device in such a way that your calves can be placed on the vibrating surfaces.
5. Switch on the device by pressing the button  (3) on the device or on the remote control.
6. Set the vibration level (1 – 10) you want by pressing the buttons + (4) and – (5) on the device or on the remote control.

The device switches off automatically after approx. 15 minutes. **Please Note:** It can be switched off manually at any time using the button  on the device or on the remote control.
7. Pull the mains adapter out of the plug socket and disconnect the barrel connector from the device.

Changing the Battery of the Remote Control

Before inserting the battery, check whether the contacts on the remote control and on the battery are clean, and clean them if necessary.

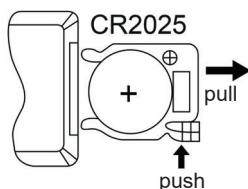1. The battery compartment is located on the narrow side of the remote control. Pull it out as shown.
2. Take out the used battery and insert a new battery (3 V, CR2025). The positive pole must be facing upwards.
3. Push the battery compartment back into the remote control to the limit stop. You should feel and hear it click into place.

Replace the battery when the power of the remote control starts to wane.

Cleaning and Storage

PLEASE NOTE!

- Do not use any caustic or abrasive cleaning agents or cleaning pads to clean the device. These may damage the surfaces.
-
- Make sure that the device is not connected to the mains power.
 - Wipe down the device with a slightly damp cloth if necessary. Then wipe it down with a dry cloth.
 - Keep the device in a dry place which is protected from sunlight and is not accessible to children or animals.

Troubleshooting

If the device does not work properly, first check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps listed, contact the customer service department.

Do not attempt to repair an electrical device yourself!

Problem	Possible cause	Solution
The device does not work.	The mains adapter is not inserted correctly in the plug socket.	Push the mains adapter into the plug socket as far as it will go.
	The device is not connected to the mains adapter correctly.	Correct the connection.
	The plug socket is defective.	Try another plug socket.
	There is no mains voltage.	Check the fuse of the electrical distributor (fuse box).
	The button  (3) has not been (correctly) pressed.	Press the button  (3) again.
The device does not react to the remote control.	The battery of the remote control is too weak.	Replace the battery of the remote control.
	There are obstructions between the device and remote control (no visual contact).	Point the remote control directly at the device.

Technical Data

Article number: 05584
ID of operating instructions: Z 05584 M DS V2 0220 dk

Device

Model: F-26
Power supply: 24 V DC
Power: max. 50 W
Vibration levels: 10

Mains adapter

Model: XH2400-2000WG
Input voltage: 100 – 240 V ~ 50 / 60 Hz, 1.5 A
Output voltage: 24 V DC, 2 A
Protection class: II

Remote control

Power supply: 3 V (1 x button cell CR2025)

Disposal

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

The battery should be removed prior to disposal of the remote control and disposed of separately from it at a suitable collection point.

This device is governed by the European Directive 2012/19/EU on waste electrical and electronic equipment (WEEE). Do not dispose of this device as normal domestic waste, but rather in an environmentally friendly manner via an officially approved waste disposal company.

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All rights reserved.
DS Produkte GmbH